

S-BASIC

AUDIOMETRIC BOOTH



Frequency 8000Hz
43.3dB



Quick and easy
assembly



Compatibility with the
main audiometers of the
market



Folding table

S-BASIC · audiometric booth

Description

The audiometric test booth **S-BASIC** uses the latest sound insulation materials. This has made it possible to develop a range of high quality booths at a very economical price, reducing their weight and obtaining a high insulation index with great functionality and aesthetic beauty, in addition to their compatibility with the main audiometers on the market.

Its exterior finish is made of washable panels, its interior looks elegant and comfortable. It has internal lighting.

Features

	S-BASIC	S-BASIC B	S-BASIC C
External booth	96x96x208 cm	110x110x208 cm	120x120x208 cm
Internal booth	86x86x190 cm	101x101x193 cm	111x111x193 cm
Window		70 x 58 cm	
Door		80 x 188 cm	
Foldable table		80 x 40 cm	
External material		Melamine	
Internal materia		Pyramid foam	
Int. illumination		Circular LED lamp	
Floor		Washable polyester	
Alimentación	220 V 50/60 Hz (Other optional voltages)		
Connection boards	5 jacks de Ø 6,3mm y 1 jack de Ø 3,5mm		
Weight	<200 Kg	<220 Kg	<240 Kg


Regulations

- The Audiometric Test Booth S-BASIC complies with:
- EN 60601-1:2006+AC:2010+A1:2013+AC:2014
- EN 60601-1-6:2010 +A1:2015
- EN 62366-1:2015 + A1:2020
- EN ISO 14971:2019 + A11:2021
- EN ISO 20417:2021
- EN ISO 15223-1:2021

This product has been validated in laboratory according to UNE-EN ISO 11957:2009 regulation and calculation of the global index of acoustic reduction, has been carried out according to the UNE EN ISO 717-1: 1997 standard UNE EN ISO 717-1:1997.

The process of designing, production and customer service will be made according to ISO 9001:2015, EN 13485:2016+AC:2018 quality standards and to the Medical device regulation (EU) 2017/745.

Noise attenuation

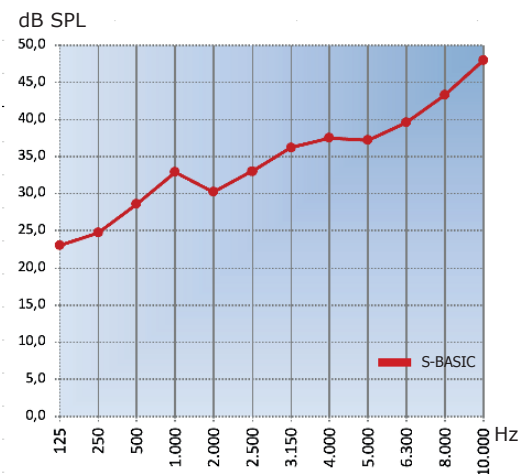


Certified test done by an Applus authorized laboratory

D'p.w:31dB

AUDIOMETER SIBEL SOUND DUO

- Two digital channels audiometer
- Color touch screen.
- Portable, with carrying bag.
- 2 models: A-AOM.
- Database: > 1000 tests with unlimited patients.
- Automatic Hughson-Westlake Audiometry.
- Data export (via Ethernet).
- USB and Ethernet connectivity.
- W50 audiometry software.

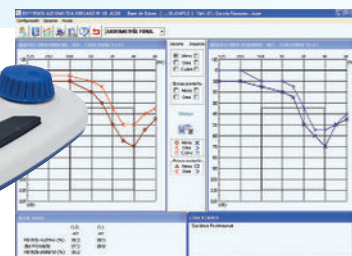

Raimon Casellas, 6-8. 08205-Sabadell (Barcelona). Spain.
National Sales: Tel. +34 93 436 00 08 e-mail: comercial@sibelmed.com
International Sales: Tel.+34 93 436 00 07 e-mail: export@sibelmed.com
Technical support: +34 93 433 54 50 e-mail: sat@sibelmed.com
www.sibelmed.com



CE



EN ISO9001:2015

SIBEL S.A.U. reserves the right to modify the technical specifications without prior notice

Sibelmed